


TMI Trace Minerals International Laboratory

good chemistry for better health

1711 Pearl Street; Suite 200
Boulder, CO 80302
USA

Phone: +1 (720) 325-2652
Facsimile: +1 (720) 325-2653

service@tracemin.com
<http://www.tracemin.com>

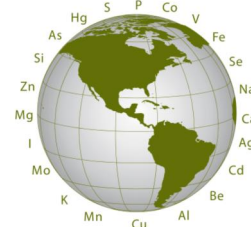

Submission Form: Hair or Nail Mineral Analysis

Requesting Clinic/Doctor:

New Customer or if contact information has changed, please fill out the fields on page 2.

Patient Name: _____

Street: _____ City: _____

State: _____ ZIP: _____ Country: _____

Phone: _____ Fax: _____

E-mail: _____

please fill out if report is to be mailed to the patient (please complete in block capitals)

Date of Birth: _____ Sex: m f

Date: _____ Patient Signature: _____

(please do not forget)

Requested Profile:

Standard Profile (P9) 35 Elements

Tested Elements:

Aluminum, Antimony, Arsenic-total, Barium, Beryllium, Bismuth, Boron, Cadmium, Calcium, Chromium, Cobalt, Copper, Germanium, Iodine, Iron, Lead, Lithium, Magnesium, Manganese, Mercury, Molybdenum, Nickel, Palladium, Platinum, Selenium, Silver, Strontium, Thallium, Tin, Titanium, Tungsten, Uranium, Vanadium, Zinc, Zirconium


Extended Profile (P10) 55 Elements

Tested Elements Parameter as profile 9 plus:

Cerium, Cesium, Dysprosium, Erbium, Europium, Gadolinium, Gallium, Iridium, Lanthanum, Lutetium, Praseodymium, Rhenium, Rhodium, Ruthenium, Samarium, Tantalum, Tellurium, Thorium, Thulium, Ytterbium

Test material:

Head hair sample (untreated) pubic hair sample beard sample Nail sample

Payment via:	Invoice to:	Doctor	Patient
Credit Card	VISA Mastercard	Card Number: _____	
valid thru (MM/YY): _____	3-digit code: _____	Signature: _____	
Bank transfer done at: _____		for \$: _____	
	Payment was made to address: service@tracemin.com		
Pre-Payment or Credit Card is Needed, otherwise samples will be held until payment is received.			

*** please turn over ***



Sampling Instructions

Hair provides a record of past and current trace element levels. As hair grows, nutrient and toxic elements are deposited from the blood stream into the hair follicle and hair shaft.

Once a trace element has been incorporated into the hair, it remains fixed. To measure these values reliably and with good reproducibility, the following criteria applies:

- Untreated head hair that is cut from the occipital area (back of head) is the preferred sample.
- Substitutes for head hair are pubic hair or clippings of finger nails (underarm hair is not acceptable).
- Because some sample loss results from thorough washing, an adequate amount of sample (about 2-3 Teaspoons filled with hair) is needed for the analysis.
- DO NOT MIX different sample types (i.e. hair with nails).
- Hair roots are NOT needed.
- Hair that has been chemically treated (“permed”, dyed, bleached, or otherwise treated) WILL NOT provide accurate results. If the hair has been treated, it should be allowed to grow for at least 2-3 months before new sampling.
- Hair analysis results are falsely elevated by:
 - Coloring or tinting agents, which often elevate a number of trace elements or toxic metals. Since today's cosmetic industry uses many different chemical agents, we have no way of accurately determining their effects.
 - Bleached and/or “Permed” hair results in altered levels of calcium, magnesium, zinc, copper, nickel and other metals.

Procedure

Hair

- Part the hair in the middle of the back of the head, and pull it out and up.
- Cut a small 1.5 to 2 inch (4.5 to 5.5 cm.) strand of hair close to the head. DISCARD ends of long strands and KEEP less than 2 inches (less than 5.5 cm) closest to the scalp.
- Place hair in sampling envelope, fill out the Sampling Instructions Form with the appropriate information and send to TMI.

Short hair

- Trim a minimum of 0.300 grams of hair from the back of the head. Use thinning scissors, if possible.
- Place hair sampling (paper) envelope, fill out the Sampling Instructions Form with the appropriate information and send to TMI.

Pubic hair

- Trim a minimum of 0.300 grams of pubic hair with normal scissors.
- Place hair in paper envelope, fill out the Sampling Instructions Form with the appropriate information and send to TMI.

Nails

- We need a minimum of 0.200 grams (200mg) of unpolished nails. This equals about ½ of a filled teaspoon. Place sample in a paper or plastic envelope, fill out the Instructions Form and send to TMI.

Send Report to:	Doctor	Patient	both addresses (\$ 5.95 surcharge)
Send Report via:	Post	E-Mail	Fax

<p>New Customer or if contact information has changed,</p> <p>Address: _____</p> <p>Phone: _____</p> <p>Fax: _____</p> <p>E-mail: _____</p>	or	<p style="text-align: center;">Clinic/Doctor Stamp</p>
Barcode H		Barcode N